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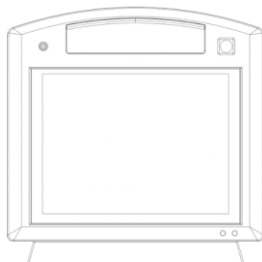
 Mui Scientific

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The Rapid Barostat Bag Pump is intended for use in the electromagnetic environment specified below. The customer or the user of the Rapid Barostat Bag Pump should assure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment – Guidance
RF Emissions CISPR 11	Group 1	The Rapid Barostat Bag Pump uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF Emissions CISPR 11	Class A	The Rapid Barostat Bag Pump is suitable for use in all establishments other than domestic, and may be used in domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes, provided the following warning is heeded: Warning: This equipment/system is intended for use by healthcare professionals only. This equipment/system may cause radio interference or may disrupt the operation of nearby equipment. It may be necessary to take mitigation measures, such as re-orienting or relocating the Rapid Barostat Bag or shielding the location.
Harmonic Emissions IEC 61000-3-2	Class A	
Voltage Fluctuations/ Flicker Emissions IEC 61000-3-3	Complies	

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The Rapid Barostat Bag Pump is intended for use in the electromagnetic environment specified below. The customer or the user of the Rapid Barostat Bag Pump should assure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	± 8 kV Contact ± 15 kV Air	± 8 kV Contact ± 15 kV Air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical Fast Transient/Burst IEC 61000-4-4	± 2 kV for power supply lines	± 2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 0.5 kV, ± 1 kV line to line & ± 0.5 kV, ± 1 kV, ± 2 kV line to ground	± 0.5 kV, ± 1 kV line to line & ± 0.5 kV, ± 1 kV, ± 2 kV line to ground	Mains power quality should be that of a typical commercial or hospital environment.
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U_T (100% dip in U_T) for 0,5 cycle 0% U_T (100% dip in U_T) for 1 cycles 70% U_T (30% dip in U_T) for 25 cycles 0% U_T (100% dip in U_T) for 5 sec	0% U_T (100% dip in U_T) for 0,5 cycle 0% U_T (100% dip in U_T) for 1 cycles 70% U_T (30% dip in U_T) for 25 cycles 0% U_T (100% dip in U_T) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Rapid Barostat Bag Pump requires continued operation during power mains interruptions, it is recommended that the Rapid Barostat Bag Pump be powered from an uninterruptible power supply or a battery.
Power Frequency Magnetic Field (50/60 Hz) IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE U_T is the a.c. mains voltage prior to application of the test level.			

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The Rapid Barostat Bag Pump is intended for use in the electromagnetic environment specified below. The customer or the user of the Rapid Barostat Bag Pump should assure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms 150 kHz to 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of the Rapid Barostat Bag Pump including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1.2\sqrt{P}$
Radiated RF IEC 61000-4-3	6 Vrms ISM/Amateur Radio bands inside 150 kHz to 80 MHz 3 V/m 80 MHz to 2,7 GHz RF communication equipment inside 80 MHz to 6 GHz	6 Vrms ISM/Amateur Radio bands inside 150 kHz to 80 MHz 3 V/m 80 MHz to 2,7 GHz RF communication equipment inside 80 MHz to 6 GHz	$d = 1.2\sqrt{P}$ 80 MHz to 800 MHz $d = 2.3\sqrt{P}$ 800 MHz to 2.7 GHz where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey ^a should be less than the compliance level in each frequency range ^b
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
a) Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Rapid Barostat Bag Pump is used exceeds the applicable RF compliance level above, the Rapid Barostat Bag Pump should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Rapid Barostat Bag Pump			
b) Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

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Recommended separation distances between portable and mobile RF communications equipment and the Rapid Barostat Bag Pump			
The Rapid Barostat Bag Pump is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Rapid Barostat Bag Pump can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Rapid Barostat Bag Pump as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1.2\sqrt{P}$	80 MHz to 800 MHz $d = 1.2\sqrt{P}$	800 MHz to 2.7 GHz $d = 2.3\sqrt{P}$
0.01	0.12	0.12	0.24
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

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Standard/Method	Description	Class/Level
CISPR 11	Power Line Conducted Emissions	Class A
CISPR 11	Radiated Emissions	Class A
IEC 60601-1-2 / IEC 61000-4-2	Electro-Static Discharge	±8kV Contact ±15kV Air
IEC 60601-1-2 / IEC 61000-4-3	Radiated Field Immunity	3 V/m, 80 MHz – 2.7 GHz & Proximity Fields of Table 9
IEC 60601-1-2/ IEC 61000-4-4	Electrical Fast Transients	±2kV – Mains
IEC 60601-1-2/ IEC 61000-4-6	Conducted RF Immunity	3 Vrms, 150 kHz – 80 MHz & 6 Vrms ISM
IEC 60601-1-2 / IEC 61000-4-39	Immunity to Proximity Magnetic Fields	134.2kHz, 2.1kHz PM, 65A/m 13.56MHz, 50kHz PM, 7.5A/m

Standard	Name
IEC 60601-1-2:2014 A1:2020	Medical Electrical Equipment – Part 1-2: General Requirements for Basic Safety and Essential Performance – Collateral Standard: Electromagnetic Compatibility
CISPR 11:2015 A2:2019	Industrial, Scientific and Medical (ISM) Radio-Frequency Equipment – Electromagnetic Disturbance Characteristics – Limits and Methods of Measurement
IEC 61000-4-2:2008	Testing and Measurement Techniques - Electrostatic Discharge Immunity Test
IEC 61000-4-3:2006 A2:2010	Testing and Measurement Techniques - Radiated, Radio-Frequency, Electromagnetic Field Immunity Test
IEC 61000-4-4:2012	Testing and Measurement Techniques - Electrical Fast Transient/Burst Immunity Test
IEC 61000-4-6:2013	Testing and Measurement Techniques - Immunity to Conducted Disturbances, Induced by Radio-Frequency Fields
IEC 61000-4-39:2017	Testing and Measurement Techniques – Radiated Fields in Close Proximity Immunity Test
ISO/IEC 17025:2017	General Requirements for the Competence of Testing and Calibration Laboratories

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IEC 60601-1-2:2014 Table 9 - Test Levels					
Test Frequency (MHz)	Band (MHz)	Service / Technology	Immunity Test Level (V/m)	Modulation ^{a)}	Test Distance
385	380 – 390	TETRA 400	27	Pulse 18 Hz	1.1m
450	430 – 470	GMRS 460, FRS 460	28	Pulse 18 Hz ^{b)}	1.1m
710	704 - 787	LTE Band 13, 17	9	Pulse 217 Hz	1.1m
745					
780					
810	800 – 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	28	Pulse 18 Hz	1.1m
870					
930					
1720	1700 - 1990	GSM 1800, CDMA 1900, GSM 1900, DECT, LTE Band 1, 3, 4, 25, UMTS	28	Pulse 217 Hz	2.2m
1845					
1970					
2450	2400 – 2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	28	Pulse 217 Hz	2.2m
5240	5100 – 5800	WLAN, 802.11 a/n	9	Pulse 217 Hz	2.2m
5500					
5785					
a) The carrier is modulated using a 50% duty cycle square wave signal. b) 18Hz pulse modulation is used as a worst-case scenario compared to FM modulation as allowed by IEC 60601-1-2:2014.					

IEC 60601-1-2 Test Levels – Immunity to Radiated Fields in Close Proximity					
Frequency	Field Strength	Modulation	Source	Environment	Test Window Size/ Distance (mm)
134.2 kHz	65 A/m	Pulse Modulation 50% Duty Cycle 2.1kHz	RFID Readers	All Environments	100x100/50
13.56 MHz	7.5 A/m	Pulse Modulation 50% Duty Cycle 50kHz	RFID Readers	All Environments	80x80/50

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